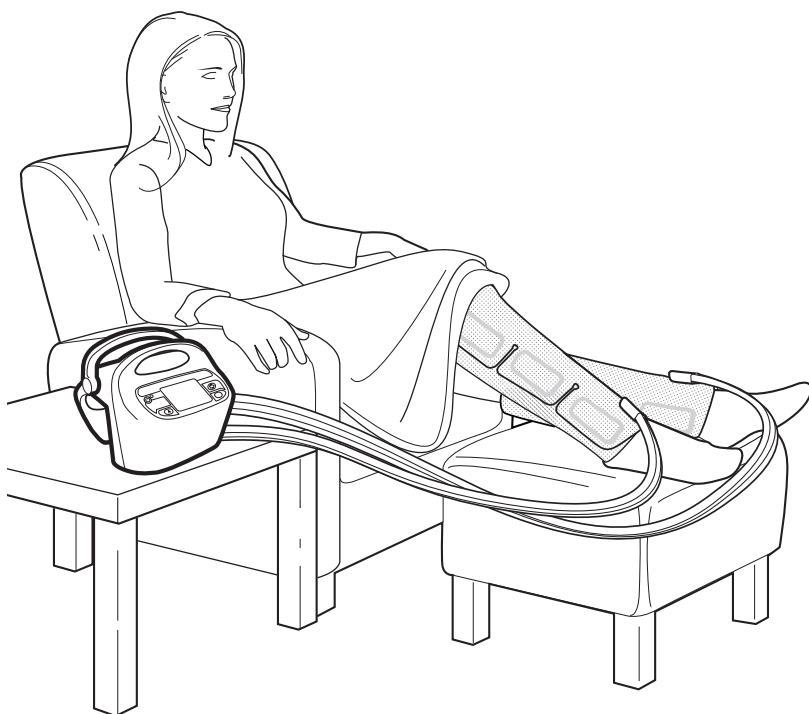
User Manual

Kendall SCD™

700 Sequential Compression System

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TABLE OF CONTENTS

Prescribing Information

Indications	1
Contraindications	2
Warnings & Cautions	3
Kendall SCD™ 700 Series Symbols Used	4

Section 1 - The Kendall SCD™ 700 Sequential Compression System

About the Kendall SCD™ 700 Series & DVT	6
What is DVT?	6
How Does the Kendall SCD™ 700 Series Help?	6
Who Should Use It?	6
Who Should NOT Use It?	6
Kendall SCD™ 700 Series Components	7

Section 2 - Using the Kendall SCD™ 700 Series

Getting Started	10
Step 1. Set Up the Controller.....	10
Step 2. Plug In	11
Step 3. Attach Tubing to Controller	12
Step 4. Apply Garment(s) to Limb(s)	13
Step 5. Starting Therapy	15
Therapy	16
During Therapy	16
Ending Therapy	17
Back-up Battery Power	18
Charging the Battery	18

Section 3 - Troubleshooting

Error Messages	20
Garment Type (Leg Sleeve or Foot Cuff)	20
Low Battery	21
Garment Connection (E12)	21
Pressure Problem (E1)	21
Pressure Problem (E2, E3, E4)	22
Pressure Problem (E8)	23
Temperature Problem (E9)	23
Call for Service (E5, E6, E7, E10, E13)	24

Section 4 - Cleaning, Contact, and Specifications

Cleaning the Kendall SCD™ 700 Series	26
Contact	27
Specifications	27

Indications

The Kendall SCD™ 700 Sequential Compression System (hereby referenced as “Kendall SCD™ 700 Series”) is designed to apply intermittent pneumatic compression to increase venous blood flow in at-risk patients in order to help prevent deep vein thrombosis and pulmonary embolism. The Kendall SCD™ 700 Series is a prescription device for use in a clinical setting or in the home. The System consists of the controller, the tubing sets (provided with the controller) and single-patient use garments (purchased separately from this controller). The garments, both leg sleeves and foot cuffs, compress the limbs to enhance venous blood movement. After the compression cycle has reached set pressure, the Controller measures the time it takes for the limbs to refill with blood and waits that period of time before the next compression is initiated.

Leg Compression

The use of the Kendall SCD™ 700 Series Compression System with Leg Sleeves is indicated for:

1. Deep vein thrombosis and pulmonary embolism prophylaxis.

Foot Compression

The use of the Kendall SCD™ 700 Series Compression System with Foot Cuffs is indicated for:

- | | |
|--------------------------------------|--|
| 1. Circulation enhancement. | 5. Extremity pain incident to trauma or surgery. |
| 2. Deep vein thrombosis prophylaxis. | 6. Leg Ulcers. |
| 3. Edema - Acute. | 7. Venous stasis / venous insufficiency. |
| 4. Edema - Chronic. | |

If you need further information regarding the Kendall SCD™ 700 Series Compression System or its clinical benefits, please contact your medical equipment provider.

Contraindications

Leg Compression

The Kendall SCD™ 700 Series Compression System may not be recommended for use with Leg Sleeve on patients with the following:

1. Any local leg condition in which the sleeves may interfere, such as: (a) dermatitis, (b) vein ligation [immediate postoperative], (c) gangrene, or (d) recent skin graft.
2. Severe arteriosclerosis or other ischemic vascular disease.
3. Massive edema of the legs or pulmonary edema from congestive heart failure.
4. Extreme deformity of the leg.
5. Suspected pre-existing deep venous thrombosis.

Foot Compression

The Kendall SCD™ 700 Series Compression System may not be recommended for use with Foot Cuffs on patients with the following:

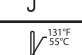
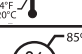
1. Conditions where an increase of fluid to the heart may be detrimental.
2. Congestive heart failure.
3. Pre-existing deep vein thrombosis, thrombophlebitis or pulmonary embolism.

Warnings & Cautions

WARNING
Do not start therapy without having completed training with your healthcare or medical equipment provider.
If you develop a skin rash or irritation near the site of the garment, stop using the system and contact your physician.
If you experience numbness, tingling or leg pain, stop using the system and contact your physician.
If you notice swelling or pain in your limbs, stop using the system and contact your physician.
Do not attempt to repair or replace broken tubing connectors as hazardous inflation of the sleeves may occur.
Explosion hazard. Not suitable for use in the presence of a flammable anesthetic mixture with air or oxygen or nitrous oxide.
Do not open battery, dispose of in fire, or short circuit. Doing so may cause the battery to ignite, explode, leak, or become hot and cause personal injury.
Do not touch or ingest any leaking electrolyte. If contact occurs, rinse skin and/or eyes immediately and seek medical attention if irritation develops. If ingested, contact local poison control center.
Use with caution on the infected or insensitive extremity.
Contact your physician if therapy is disrupted for an extended period of time.
If you experience pain or moderate discomfort associated with the therapy, please remove garments, discontinue use and contact your physician immediately.

CAUTION
Federal (U.S.A.) law restricts this device to sale by or on the order of a physician.
Patients with diabetes or vascular disease require frequent skin assessment.
Do not operate the Controller if the power cord is damaged.
Do not operate near water or if wet.
Do not operate while walking.
Do not operate while driving.
Keep Kendall SCD™ 700 Series system away from children and pets.
Contact your medical equipment provider if you are unable to operate the system.

Kendall SCD™ 700 Series Symbols Used

	Federal (USA) law restricts this device to sale by or on the order of a physician.
	Caution consult accompanying documents
	Consult instructions for use
	Non Sterile
	Catalog number
	Type BF protection against electrical shock
	Manufacturing date code
	Manufacturer
	Latex-free
	Serial Number
	Class II Equipment
	Protection against fluid ingress; spraying water
	WEEE (Waste from electrical and electronic equipment)
	Keep away from sunlight
	Keep Dry
	Store in temperatures between -4°F and 131°F (-20°C and 55°C) Operate in temperatures between 50°F and 104°F (10°C and 40°C)
	Humidity limitations
	Classified with respect to electrical shock, fire and mechanical hazards only in accordance with UL60601-1 CAN/CSA C22.2 NO.601.1

Section 1

The Kendall SCD™ 700 Series Sequential Compression System

About the Kendall SCD™ 700 Series & DVT

What is DVT?

Deep Vein Thrombosis or DVT is the formation of a blood clot in a deep vein in the legs. Blood clots may form in the deep veins of the legs following surgery or periods of time when you are unable to walk. If untreated, blood clots can cause disability or even become life threatening. One of the most serious complications that can happen with DVT is that a blood clot could travel to the lungs, block blood flow and cause pressure on the heart muscle (known as a pulmonary embolism).

How Does the Kendall SCD™ 700 Series Help?

The Kendall SCD™ 700 Series Controller compresses air at regular intervals into leg and foot garments wrapped around your legs or feet. The air gently squeezes your limbs, one at a time, beginning at the bottom of the garment. Therapy is customized for you based on the time it takes your veins to refill after squeezing. After each set of compressions there is a small wait, then the compressions begin again. This pattern of compression and release increases blood flow to the legs, helping to prevent DVT and pulmonary embolism.

Who Should Use It?

You should ONLY use the Kendall SCD™ 700 Series if it is prescribed for you by your physician.

Your physician will determine if you need this therapy, and if the therapy is safe for you to use based on your health and any existing conditions that you may have.

Review your prescription prior to use. If you have any questions, contact your physician.

Who Should NOT Use It?

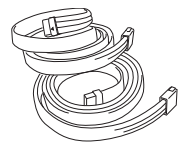
If you have certain health conditions, this therapy may be unsafe for you. Do NOT use the Kendall SCD™ 700 Series unless it has been prescribed for you by your physician.

Kendall SCD™ 700 Series Components


Components

Controller & Tubing Set

Controller

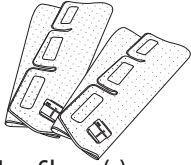
Tubing

Garments (packaged separately)

Leg Sleeve(s)

or

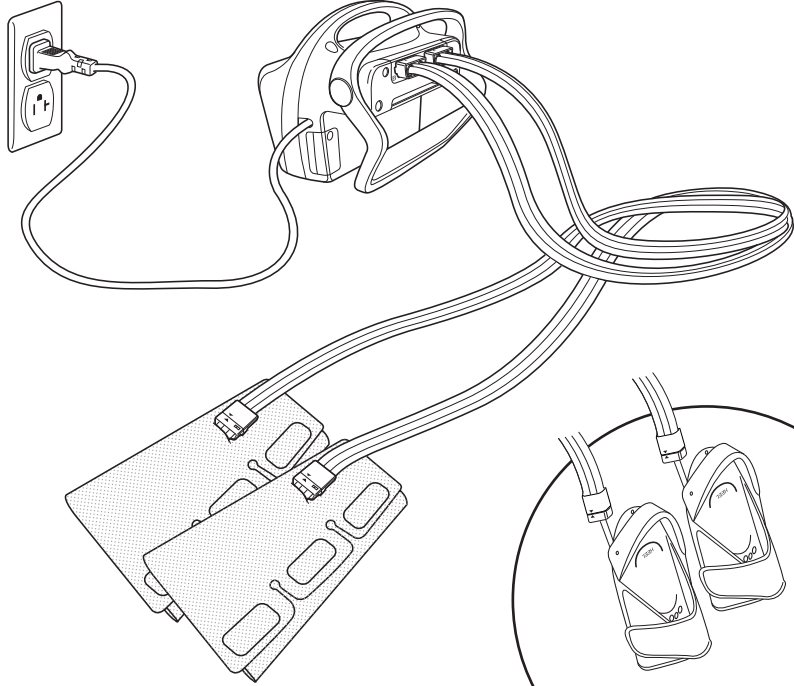
Foot Cuff(s)

NOTE

Inspect for missing, soiled, or damaged components prior to use. For replacement components, contact your medical equipment provider.

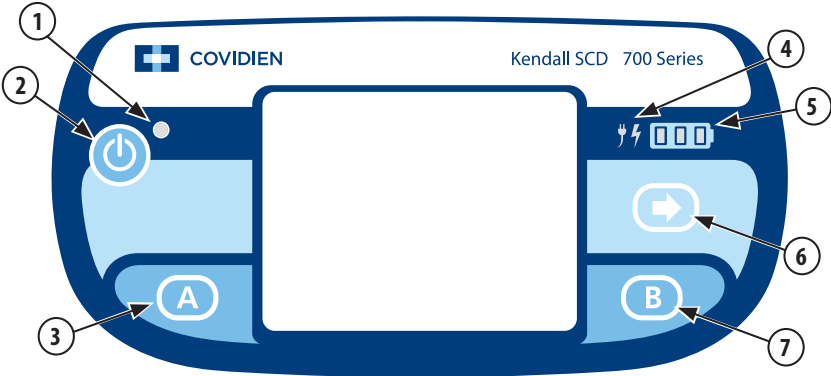Controller Front Panel Display

Item	Explanation	Item	Explanation
1	Power On Indicator	5	Battery Status Indicator
2	Power On/Standby Button	6	Arrow Button (For Service Personnel Only)
3	A - Button	7	B - Button
4	AC Power/Battery Charging Indicator		

Section 2

Using the Kendall SCD™ 700 Series

Getting Started

Step 1. Set Up the Controller

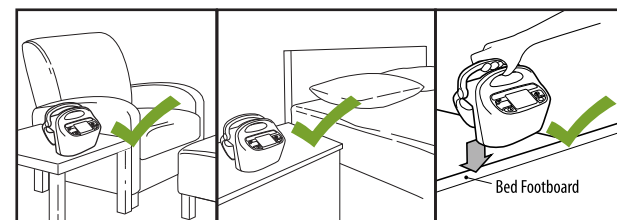
Place the Controller near where you will be sitting or resting and within reach of an AC power outlet. You can set the Controller on a bedside table or other flat surface, or you can clip the Controller onto a foot board.

Placing on a Flat Surface

- Place the Controller on a flat, horizontal surface, such as a table, reasonably close to where you will be sitting or resting.
- Make sure the Controller is not resting on any cords or tubing.
- Do not cover the Controller with blankets or other materials that may block the air flow around the unit.

To Attach to a Foot Board

- To open the gap, using one or two hands, grasp and squeeze as shown in the picture.
- Lower the Controller onto the foot board so it straddles the foot board and release the bed clamp.
- Make sure the Controller is firmly fastened.

NOTE

Save packaging materials for return of Controller and tubing after therapy.

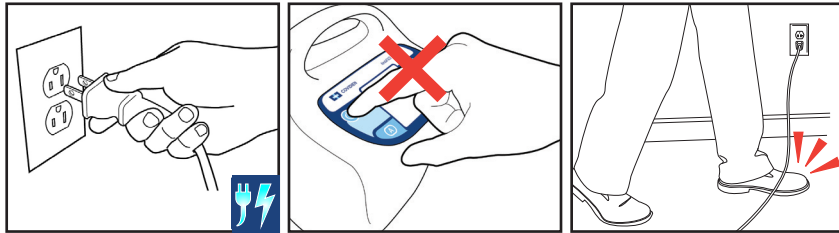
Do not cover Controller with blankets or other materials that may block the air flow around the unit.

WARNING

To avoid danger of suffocation, keep packaging material away from children.

Step 2. Plug In

- Plug the Controller power cord into an AC power outlet. The blue AC Power Indicator will illuminate.
- Do NOT turn on the Controller.

NOTE

If you turn the Controller on before completing steps 3 and 4 you will see a garment connection error. You will have to turn off the Controller and complete steps 3 and 4 to correct the error. Refer to Section 3, Error Screens.

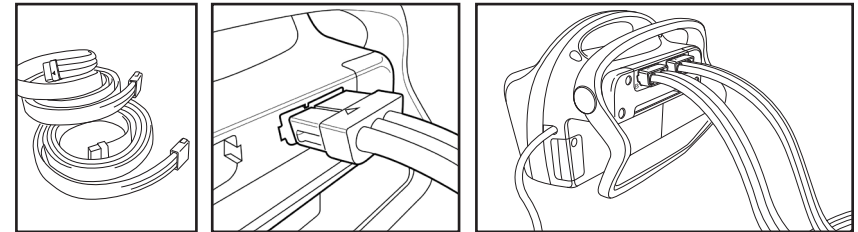
CAUTION

To avoid tripping hazard, do not place power cord, tubing, or Controller in walk areas.

Step 3. Attach Tubing to Controller

The Controller can operate with one or two garments attached.

- Tubing should be fully unwound before connecting to Controller and garments.
- Plug the tubing set(s) into the connection point(s) on the back of the Controller.
- Route the tubing from the Controller toward your limb(s). To avoid tripping hazard, do not place power cord, tubing, or Controller in walk areas.

NOTE

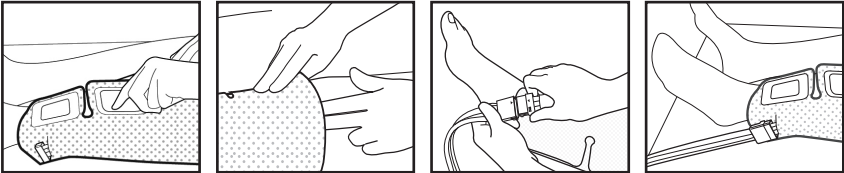
Match the port on the left, marked A, with the left limb. Match the port on the right, marked B, with the right limb. Although not required, this matching can make it easier to troubleshoot any problems that may occur.

Step 4. Apply Garment(s) to Limb(s)


The Controller can operate with one or two garments attached. Consult the instructions provided with your garment(s) to make sure each garment is aligned correctly on your limb(s).

Applying Leg Sleeve(s)


- Open your sleeve with tubing side down (away from your limb) and place under your leg.
- Note that the tubing connector should be at the ankle.
- Wrap the garment around the limb and press to secure.
- Check the fit to be sure the garment is not too loose or too tight. You should be able to fit two fingers between the garment and your limb.
- Plug the other end of the tubing into the connection point(s) on your garment(s).
- Make sure the tubing set(s) are securely connected at both ends and there are no kinks in the tubing.

Thigh length:

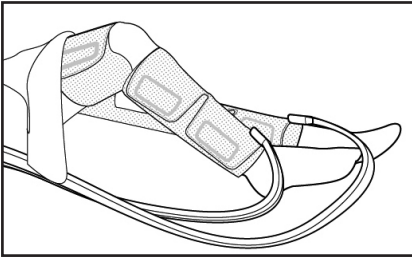
Knee length:

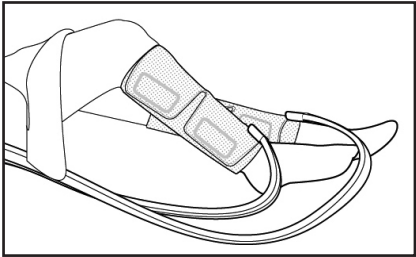WARNING
Make sure garments are positioned so that the tubing does not form pressure points on the limb(s).
Remove garment if you experience numbness, tingling, or leg pain and call your physician.
The tubing connector should be at the ankle.

Applying Foot Cuff(s)

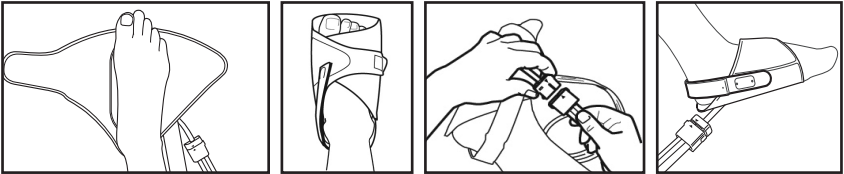- Open your garment with tubing side down (away from your limb) and place under your foot.
- For proper fit, ensure the heel is aligned with the edge of the cuff.
- Use the image on the cuff for help aligning.
- Wrap the garment around the foot and hook securely.
- Wrap the strap around the back of your ankle, attaching it to the front of the cuff.
- The cuff should fit snug, but not too tight.
- Plug the other end of the tubing into the connection point(s) on your garment(s).
- Make sure the tubing set(s) are securely connected at both ends and there are no kinks in the tubing.

WARNING
Make sure garments are positioned so that the tubing does not form pressure points on the limb(s).
Remove garment if you experience numbness, tingling, or leg pain and call your physician.
The tubing connector should be at the ankle.

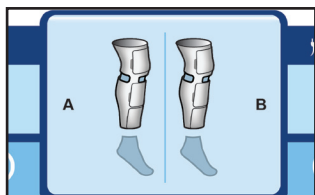
Step 5. Starting Therapy

- Press the Power On/Standby button to turn on the Controller. Images on the display will blink while the Controller completes the Garment Detection procedure and internal checks.
- Wait for the Controller to complete the internal checks.


If You Are Using Leg Sleeve(s)

- If using leg sleeves, the control panel will show an image of the limb(s) receiving therapy.
- When the leg sleeve image(s) appear, leg compression therapy has begun.

If You Are Using Foot Cuff(s)

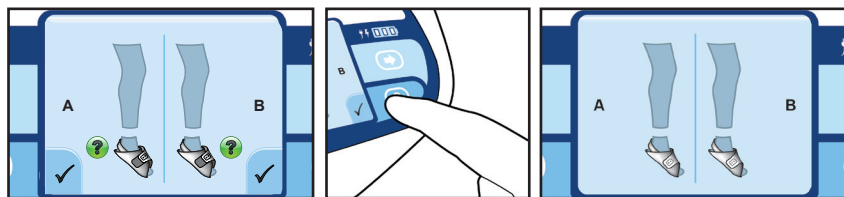- If using foot cuff(s), the Controller will beep, displaying a check mark next to the A and B button(s) corresponding to the foot wrapped with a foot cuff garment.
- When the check mark(s) appear, press the A and/or B button(s) to confirm that you are using a foot cuff on that limb.
- When the check mark(s) disappear and the controller stops beeping, foot compression therapy has begun.

Therapy

During Therapy

Treatment

Follow your physician's instructions for therapy.

Check Skin Regularly

Check the skin under the leg sleeve(s) or foot cuff(s) at least every 8 hours. Look for any changes in skin conditions. If you see any changes in skin conditions, remove the garments and call your physician.

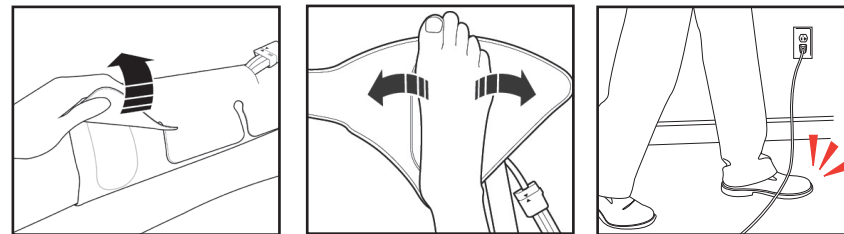
NOTE

You can wear the leg sleeve(s) or foot cuff(s) on bare skin or over light clothing or stockings.

Remove Before Walking or Bathing

Remove the leg sleeve(s) or foot cuff(s) when you need to get up and walk or bathe:

- Press the Power On/Standby button to turn off the Controller.
- Unwrap the garment(s) from your limb(s),
- Leave the tubing attached to garments and Controller.

CAUTION

To avoid tripping hazard, do not place power cord, tubing, or Controller in walk areas.

Do not attempt to walk while wearing garment(s).

Be cautious of loose, dangling tubing or cord when attempting to walk. Do not wrap tubing or cord around neck.

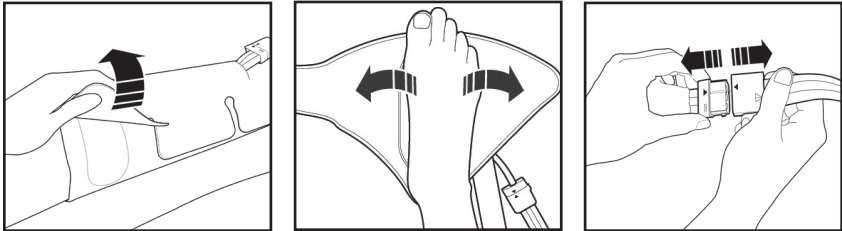
To Resume Therapy

- Reapply leg sleeve/foot cuff(s). Remember to keep the tubing free of kinks. To avoid tripping hazard, do not place power cord, tubing, or Controller in walk areas. Refer to Section 2, Getting Started, Steps 4 and 5.
- Press the Power On/Standby button to turn on the Controller.
- Remember if you are using Foot Cuffs, you will have to press the indicated A and/or B button(s) to confirm.

Ending Therapy

When you have completed your therapy according to your physician’s directions:

- Press the Power On/Standby button to turn off the Controller.
- Unwrap and remove the garment(s) from your limb(s).
- To disconnect the tubing from garment(s), grasp opposing sides of the tubing connector assembly firmly in each hand and pull apart.
- Discard garment(s).
- To disconnect the tubing from the Controller, grasp the Controller firmly in one hand and the connector in the other hand and pull apart.
- Remove the power cord from the electrical outlet by grasping the connector. Do not pull the cord.
- Return the system to your medical equipment provider. Follow provider’s directions for repackaging Controller and tubing for return.

NOTE

Do not dispose of Controller or tubing. Return Controller and tubing to your medical equipment provider.

If something is spilled, or the Controller and tubing become dirty, refer to Section 4 - Cleaning.

WARNING

Do not walk or bathe while using device.

Back-up Battery Power

It is recommended that you connect the Controller to AC power, while on or on standby. If you are unable to use AC power, the Controller can run on battery power.

Fully charged

15 - 40 minutes of charge

Less than 15 minutes of charge

When there is less than 15 minutes of battery charge remaining, the Controller will beep and the low battery screen will display.

Charging the Battery

The battery will begin charging as soon as the Controller is plugged into an AC power source.

Battery charging

The battery pack contains Lithium Ion (Li-Ion) battery cells and must be used properly for safety and to maintain optimal performance.

WARNING

Do not open battery, dispose of in fire, or short circuit. Doing so may cause the battery to ignite, explode, leak, or become hot and cause personal injury.

Do not touch or ingest any leaking electrolyte. If contact occurs, rinse skin and/or eyes immediately and seek medical attention if irritation develops. If ingested, contact local poison control center.

CAUTION

Store Controller between -4°F and 131°F (-20°C and 55°C).

Do not drop, impact, or immerse in water.

In the event of a battery problem, call your medical equipment provider.

Section 3

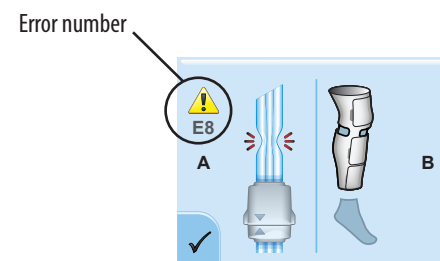
Troubleshooting

Error Messages

This section describes the errors that may occur and what you can do to troubleshoot any problems.

Note that your screen may look slightly different, depending on whether you are using one or two leg sleeves or foot cuffs.

To find an error number, look in the upper left corner of the screen (not all screens have an error number).

NOTE

If you are unable to resolve a problem, turn off the device and call your medical equipment provider.

Garment Type

Error Screen	Description	Troubleshooting
	<p>The garment type does not match the user-selected type (foot cuff or leg sleeve).</p>	<p> Press the A and/or B button(s) next to the check mark(s) to change the garment selection.</p>

If you are unable to resolve a problem, turn off the device and call your medical equipment provider.

Low Battery

Error Screen	Description	Troubleshooting
	Less than 15 minutes of battery charge remaining.	Plug the Controller into an AC power outlet.

Garment Connection (E12)

Error Screen	Description	Troubleshooting
	No leg sleeve(s) or foot cuff(s) are detected.	<ul style="list-style-type: none"> Check for disconnected leg sleeve(s) or foot cuff(s) or tube sets and reconnect. Press the A and/or B button(s) next to the check mark

Pressure Problem (E1)

Error Screen	Description	Troubleshooting
	The pressure is too high.	<ul style="list-style-type: none"> Check that the tubing is not blocked or kinked. Check that the leg sleeve(s) or foot cuff(s) are not kinked or pinched. Press the Power button to turn Controller to standby mode, and then press again to restart.

If you are unable to resolve a problem, turn off the device and call your medical equipment provider.


Pressure Problem (E2, E3, E4) continued

 	<p>The pressure is not correct.</p>	<ul style="list-style-type: none"> Check for disconnected leg sleeve(s) or foot cuff(s) and reconnect. Check that the tubing is not blocked or kinked. Check for leaks in the sleeve or the tube connections. Check that the leg sleeve(s) or foot cuff(s) are not kinked or pinched. After addressing any problems, press the A and/or B button(s) next to the check mark(s).
--------------	-------------------------------------	--


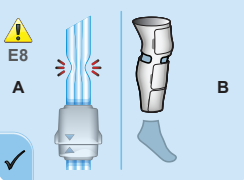
If you are unable to resolve a problem, turn off the device and call your medical equipment provider.

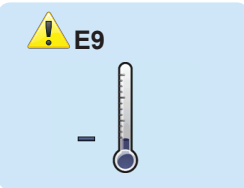Pressure Problem (E8)

	<p>The pressure is too high.</p>	<ul style="list-style-type: none"> Check that the tubing is not blocked or kinked.  Press the A and/or B button(s) next to the check mark
--	----------------------------------	--

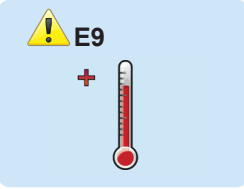Temperature Problem (E9)

Error Screen	Description	Troubleshooting
	<p>The temperature of the Controller is too low.</p>	<ul style="list-style-type: none">  Press the Power button to turn Controller to standby mode. Allow the Controller to warm to room temperature.  Press the Power button to turn the Controller back on.
	<p>The temperature of the Controller is too high.</p>	<ul style="list-style-type: none">  Press the Power button to turn Controller to standby mode. Make sure the Controller is not covered by bedding and that the fan port, located near the power cord is not obstructed. Allow the Controller to cool to room temperature.  Press the Power button to turn the Controller back on.

If you are unable to resolve a problem, turn off the device and call your medical equipment provider.

“Call for Service” Error (E5, E6, E7, E10, E13)

Error Screen	Description	Troubleshooting
	<p>The Controller has an error that requires service.</p>	<ul style="list-style-type: none">  Press the Power button to turn Controller to standby mode.
		<ul style="list-style-type: none">  Call your medical equipment provider.
		
		
		

Section 4

Cleaning, Contact, and Specifications

Cleaning the Kendall SCD™ 700 Series

CAUTION

Make sure that the Controller is in standby mode and that the power cord is disconnected from an outlet before attempting to clean.

Do not attempt to open the Controller.

Do not immerse the Controller in water or any liquid.

Do not immerse the tubing in water or any liquid.

Avoid excessive spraying, especially in the areas of the connection ports on the back of the device. If any liquid enters the ports, internal component damage will likely result.

Cleaning Controller

The outer surface of the Controller and tubing can be cleaned with a soft cloth dampened with the cleaners listed below.

Kendall SCD™ 700 Series Controller CLEANERS

0.5% bleach solution (One part bleach to ten parts water, 3-6% solution)

70% Isopropanol alcohol

3 - 7% Hydrogen Peroxide

Ensure that the connectors are fully dry before re-connecting to the Controller.

If needed, contact your medical equipment provider for replacement Controller and tubing.

Cleaning Leg Sleeves and Foot Cuffs

Do not wash leg sleeves or foot cuffs. If needed, contact your medical equipment provider for replacement garments.

Contact

The Kendall SCD™ 700 Series Controller contains no user-serviceable parts. For questions, including maintenance, please contact your medical equipment provider.

Specifications

Safety Standards	Built to UL60601-1, CSA-C22.2 No. 601.1-M90, EN60601-1, and IEC 60601-1-2 Standards UL Classified File # E189131
Device Classification	Class II Equipment Internally Powered, Portable Type BF Applied Parts Not AP or APG Equipments
Mode of Operation	Continuous
Ingress of Water Protection	IPX3 (IEC 60529)
Power Requirements	100-120 V~, 60Hz, 50VA

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**For any questions please contact your local
Kendall SCD™ 700 Series provider.**

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**Rx
ONLY**

Caution, consult
accompanying
documents

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IFU No. SMF1111-009

COVIDIEN, COVIDIEN with logo, and Covidien logo are U.S.
and internationally registered trademarks of Covidien AG.

™*Trademark of its respective owner. Other brands are
trademarks of a Covidien company.

U.S. Patents 5,876,359; D621,471S.

Covered by U.S. Patents when used with an inflatable
chamber: 6,231,532; 6,468,237; 7,076,993; 7,354,410;
7,354,411, 7,490,620.

©2011 Covidien.

Printed in China.

COVIDIEN™

 Covidien LLC, 15 Hampshire Street,
Mansfield, MA 02048 USA.

www.covidien.com

REV 07/2012